


Sure Quality is on
Sure-endo

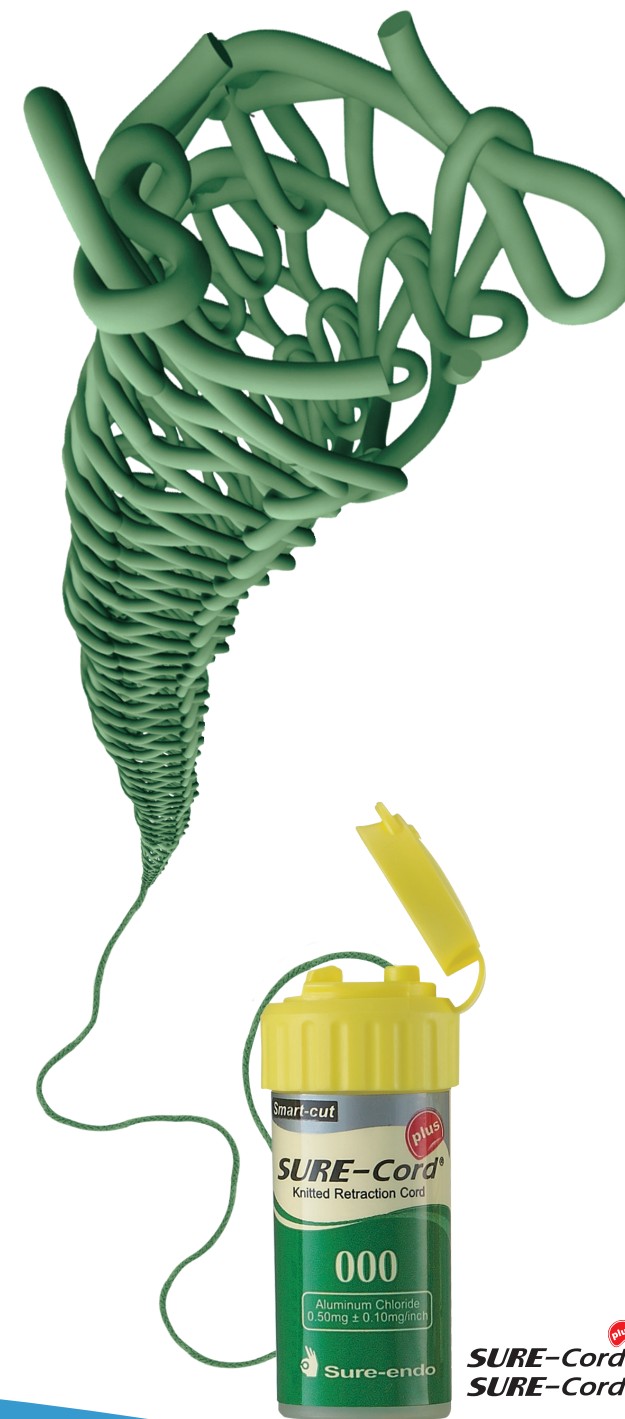

SURE-Cord^{plus}[®]

Knitted Retraction Cord

Retraction&Hemostasis cord from
Sure-endo

Just try it
feel the differences

Clean | Strong | Rapid



Sure Dent Corporation


#809, 52, Sagimakgol-ro, Jungwon-gu, Seongnam-si, Gyeonggi-do, Rep. of Korea
Tel. +82-31-737-4720 E-mail. suredent@suredent.com Web. www.suredent.com



SD-02 | Rev. 10 / 2019

SURE-Cord®

Knitted type Retraction cord

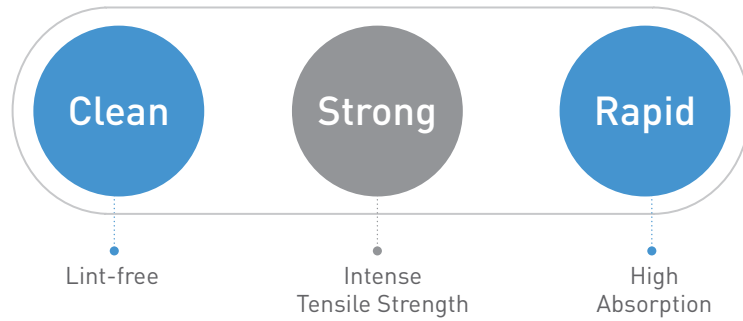

SURE-Cord® is knitted type of non-impregnated retraction cord which is made of 100% cotton.

SURE-Cord®'s knitted interlocking chains carry significantly greater quantities of hemostatic solution and provide an excellent medium for gel absorption and placement.


Advantages

- Quick tissue displacement with detailed margins and perfect impressions
- High absorption of bleeding by fine and tight knit
- Clean surface allows the cervical resin clearly
- One touch smart cutting cap ensures the clean and easy cut
- Easy placement in the sulcus
- Length indicator helps to cut the desired length
- Translucent bottle allows the easy confirmation of the remained cord

Sizes & Length


Length	100 Inch (254cm)
	#000 Green
	#00 Brown
	#0 Purple
	#1 Blue
	#2 Orange
	#3 Yellow

Application

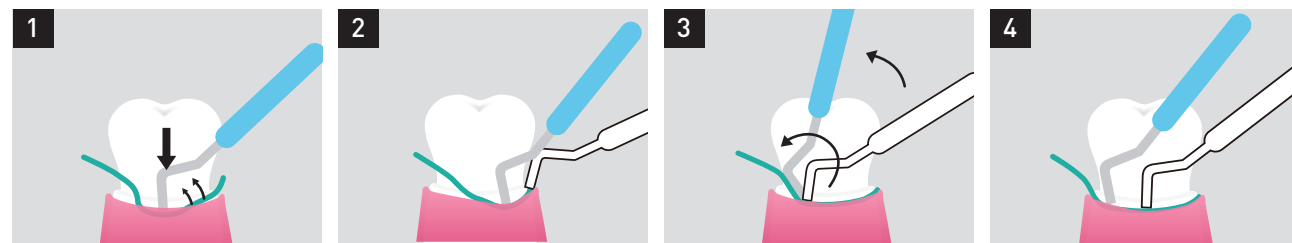

SURE-Cord® plus

Retraction & Hemostasis cord

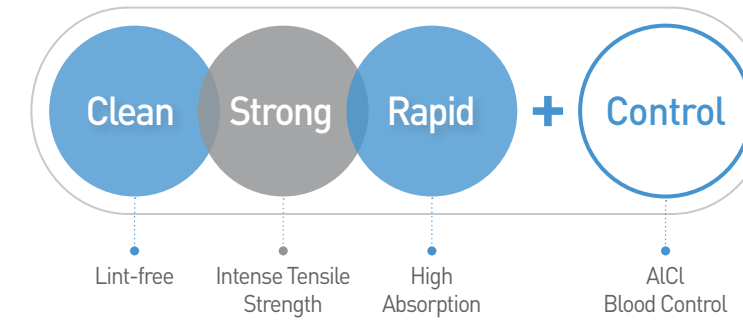

SURE-Cord® Plus is Knitted type of Retraction Cord made of Split Microfiber which is more absorbent than cotton.

SURE-Cord® Plus is impregnated with Aluminum chloride which reduces the liquid in gingival pocket and closes the smaller blood vessels.

Advantages

Impregnated with Aluminum Chloride

- Gentle astringent effect and hemostasis
- No cardiovascular problem(Epinephrine Free)
- No blackening of tissue

Intense tensile strength even with AICI

One touch cutting cap with Zirconia Blade

- Zirconia blade ensures the clean and hygienic condition without corrosion

Lint-free

- No lint provides the cleanest impression and cervical resin

Thinnest size (#000.0)

- Best application for the double cord technique


Knitted with Split Microfiber

- Perfect absorption action
- x500 magnified SEM photo

Sizes & Length

Length	100 Inch (254cm)
	#000.0 Black
	#000 Green
	#00 Brown
	#0 Purple
	#1 Blue
	#2 Orange

Knitted with Split Microfiber

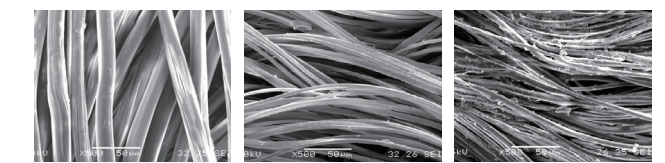

Before Split Split Microfiber After AICI Finished Products

Consultants' comments

The 000.0 size is an excellent addition to the traditional cord sizes

Sure-Cord+ is far more resistant to shredding and wrapping around a high speed bur

I love using the double- and triple-zero for two-cord technique



DENTAL ADVISOR

